Infinity

H3103

- EN USER MANUAL
- DE BEDIENUNGSANLEITUNG
- FR MODE D'EMPLOI
- NL GEBRUIKSAANWIJZING
- ES MANUAL DEL USUARIO
- PT MANUAL DO UTILIZADOR
- GR ΕΓΧΕΙΡΙΔΙΟ ΧΡΗΣΤΗ
- TR KULLANIM KILAVUZU
- IT MANUALE PER L'UTENTE
- NO BRUKSANVISNING
- SV ANVÄNDARHANDBOK
- DA BRUGERMANUAL

Silk'n Infinity

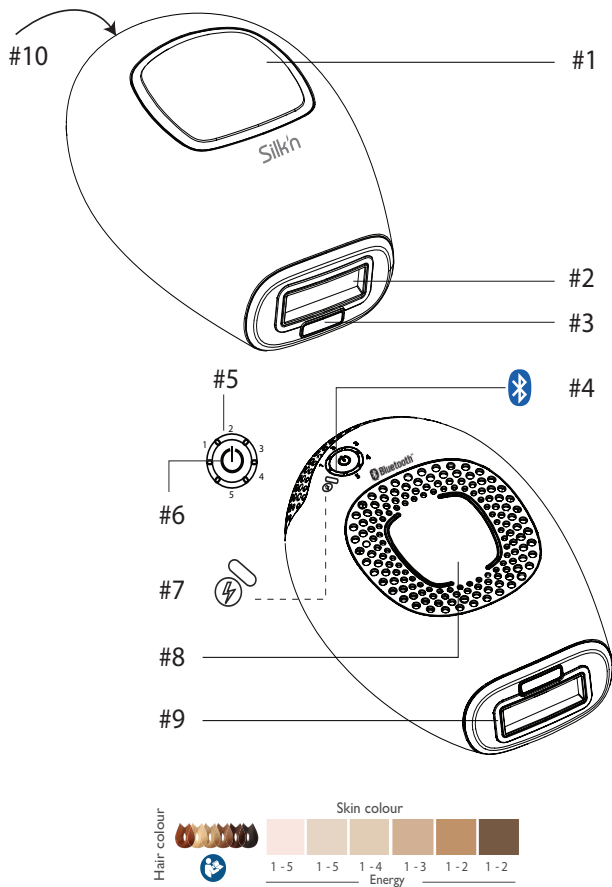
Table of contents

1. Warnings & contraindications.....	2
2. Signs and symbols	9
3. Safety with Silk'n Infinity	10
4. Get to know your Silk'n Infinity	11
5. Device controls and settings.....	12
6. What to expect with Silk'n Infinity	15
7. Possible side effects.....	15
8. Treatment plan.....	16
9. Get started!.....	18
10. Cleaning & maintenance	20
11. Connect your Silk'n Infinity to Bluetooth	20
12. Troubleshooting	21
13. Customer service.....	22
14. Warranty	22
15. Disposal	22
16. Specifications	23
17. Legal notice	24

1. Warnings & contraindications

EN

A statement or event that indicates a possible or impending danger, problem, or other unpleasant situation when using the device.

Electricity & safety warnings

The user manual is based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws. Please read all warnings and contraindications before use.

- ⦿ Always ensure that your hands are dry before connecting the adapter to the electricity supply.
- ⦿ Do not use the device while bathing or taking a shower.
- ⦿ Keep the device away from water! Do not place or store the device where it can fall, be pushed or placed into a tub, sink or any other vessel containing water.
- ⦿ Never immerse the device or any of its parts in water.

- ⦿ If the device has come into contact with water, disconnect it immediately from the electricity supply before touching it.
- ⦿ Never use the adapter if it has a damaged cable or plug, or isn't working properly, or has been submerged in water. If that happens, unplug the adapter immediately before touching it. Failure to unplug the adapter may result in an electric shock!
- ⦿ Do not operate the device in environments where there is a risk of explosion. Keep the device and all its parts away from open flames and hot surfaces.
- ⦿ Never leave the device unattended when connected to the electricity supply. Keep the adapter accessible at all times, so it can be removed from the electricity supply easily in case of emergency.
- ⦿ Always switch the device off and disconnect it from the electricity supply when you are not using the device, when you intend to clean it, or in the event of a fault.
- ⦿ Never pull the adapter out of the socket by the cable. Always remove it by taking hold of the adapter.
- ⦿ Place the adapter cable in such a way that it does not pose a tripping hazard.

- Do not kink the adapter cable and do not lay it over sharp edges.
- Keep the device away from the eyes. Possible hazardous optical radiation is emitted from this device. Do not look at the operating lamp. Eye injury may result.
- Changes or modifications to this device not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
- Do not use if the parts or accessories are not supplied by the manufacturer as defined in this manual, or if the device is damaged in any way, does not appear to be working properly, or if you see or smell smoke. In such cases, stop using the device and contact Customer Service.
- Only authorized personnel are permitted to perform repairs. Do not attempt to modify, open or repair the device. You may expose yourself to dangerous electrical components and may cause serious harm to yourself. It will also void the warranty.
- The manufacturer accepts no liability for damage or injuries caused by improper or incorrect use.
- Children must not play with the device or any of its parts.

- ⦿ Children must not play with the packaging. They may injure themselves with the packaging material, or swallow small parts and choke.
- ⦿ Do not dispose with household waste! The law requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices. For details see section 15. Disposal.

Contraindications

- ⦿ Silk'n Infinity must not be used by the following people: children, adults with limited physical, sensory or psychological capacities, or those lacking experience or knowledge in how to use the device.
- ⦿ Silk'n Infinity is not designed for everyone. Hair removal by lasers or intense pulsed light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern and South Asian heritage treated on the face and neck.
- ⦿ Silk'n Infinity is not suitable for use on children or animals.

Do not use if you:

- ⦿ are pregnant or nursing.
- ⦿ have a history of keloidal scar formation.
- ⦿ have suffered or suffer from epilepsy.
- ⦿ have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- ⦿ have a history of skin cancer or areas of potential skin malignancies.
- ⦿ have received radiation therapy or chemotherapy treatments within the past 3 months.
- ⦿ have received drugs from your physician for a skin condition within the past 6 months.
- ⦿ have any other condition which in your physician's opinion would make it unsafe for you to use the device.
- ⦿ cannot sense heat properly because of an illness, or other health related issue.

Do not use on these parts of the body:

- ⦿ On tanned skin. Refer to section 5.5 Exposure to the sun.
- ⦿ Around or near the eyes. Use below the cheekbones.

- On tattoos or permanent makeup, dark brown or black spots, (such as large freckles, moles, or warts), nipples, genitals, or lips.
- Areas with injected fillers or toxins.
- A part of the body that has a history of herpes outbreaks.
- Skin that has active eczema, psoriasis, lesions, open wounds or infections (cold sores), abnormal skin conditions caused by systemic or metabolic diseases (diabetes for example).
- Areas where you may want hair to grow later. The results are permanent.

Additional warnings

- Remove all jewellery before use.
- Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol), or acetone to clean the skin before use.
- Select the energy setting carefully. When the energy intensity is increased, hair removal improves, but the risk of side effects and complications increases.
- Follow the hair and skin colour chart when selecting energy levels to avoid the risk of burns, blisters or changes in skin colour (see the chart at the front of the manual).

- Always do a spot test before you treat a specific part of your body.
- Do not treat the same area of skin more than once per hair removal session.
- Avoid overlapping pulses or gaps when you treat an area.
- Use a lower energy level if you experience mild discomfort.
- Do not expose treated areas to the sun, and do not use energy levels 2-5 on tanned skin or after recent sun exposure (including artificial tanning machines). Refer to section 5.5 Exposure to the sun.
- Stop use immediately if you see full redness of the skin, blistering or burns.
- Stop use immediately if you have any pain.
- Stop use immediately if the skin or the device becomes too hot.

2. Signs and symbols

	<p>Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.</p> <p>Simplified Declaration of Conformity: In accordance with article 10(9) of the EU declaration of conformity, the radio equipment type Bluetooth LE is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: https://silkn.eu/media/downloads/EU-Declaration-of-Conformity-Infinity-Fast.pdf</p>
	<p>Class II. The adapter is protected against electric shock by means of double insulation.</p>
<p>IP30</p>	<p>The device is protected against the ingress of objects larger than 2.5mm diameter, such as fingers, to protect against electric shock.</p>
	<p>Do not use in wet environments.</p>
	<p>Do not dispose with household waste! The WEEE directive requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices.</p>
	<p>Conforms to technical regulations required by the Eurasian Customs Union.</p>

-  This device is not suitable for use in the presence of flammable anaesthetic mixture with air or with Oxygen or Nitrous Oxide.

3. Safety with Silk'n Infinity

EN

With Silk'n Infinity safety comes first. Enhanced Home Pulsed Light™ (eHPL) technology provides superior safety with lower energy levels.

eHPL™ is able to achieve long-term hair removal results at a fraction of the energy level used in other light-based hair removal equipment. The technology combines 2 energy sources: galvanic energy and optical energy. The galvanic energy enables better penetration of the optical energy to the hair shaft by temporarily widening the skin pores. The hair shaft absorbs the light, deactivates the hair follicle and hair growth is disabled.

The low energy used reduces its potential to cause harm or complications and contributes to your overall safety.

Silk'n Infinity protects your skin

Light-based hair removal is not suitable for tanned skin. Silk'n Infinity comes with a built-in skin colour sensor that is designed to measure the complexion of the applied surface and enables energy levels according to skin complexion.

Furthermore, the treatment surface, through which pulses of light are delivered, is recessed inside the device. This enables Silk'n Infinity to protect your skin by avoiding direct contact between the treatment surface and the skin.

Silk'n Infinity protects your eyes

Silk'n Infinity has a built-in safety feature for eye protection. It has been designed so that a light pulse cannot be emitted when the treatment surface is facing open air. The safety feature enables treatment only when the treatment surface is in contact with the skin.

4. Get to know your Silk'n Infinity

Silk'n Infinity is a powerful electrical device and should be used with special attention to safety. **Before use, it is important that you read and understand the instructions on how to use the device including after-use procedures, and follow them strictly.** We recommend you familiarise yourself with the instructions in this manual before each treatment session.

4.1 Package contents

Package contents vary according to the option purchased.

Infinity items

- Silk'n Infinity device
- Adapter
- User manual
- Warranty card
- Beauty case

4.2 Device description

Silk'n Infinity is a home-use device using eHPL™ technology for the removal of unwanted hair. The device has a treatment surface and a pulse button. It also has a pair of electrodes, base and touch, for creating and utilizing a galvanic treatment circulation. Refer to the diagram at the front of the manual to see the location of each feature.

1. Pulse button
2. Treatment surface
3. Skin colour sensor
4. Bluetooth indicator light
5. Energy level indicator lights
6. Control button
7. Ready/Skin tone warning indicator light
8. Base electrode
9. Touch electrode
10. Adapter socket

4.3 Intended use

Silk'n Infinity is intended for removal of unwanted hair. The device is also intended to effect long term or permanent hair reduction. Permanent hair reduction is defined as a long-term stable reduction in the number of hairs regrowing after a treatment regime. You may use the device to remove unwanted body hair. Ideal body areas include the underarms, bikini line, arms, legs, face (below the cheekbones), back, shoulders and chest.

Silk'n Infinity is not a medical appliance. It is intended exclusively for cosmetic treatments. The device is only intended for private use and is not suitable for commercial or therapeutic use.

Do not use for any purpose other than as described in this manual.

5. Device controls and settings

5.1 Device controls

Refer to the diagram at the front of the manual to see the location of each control.

Control	Action	Function / Mode
Control button (#6)	Short press	Each press increases the energy level. The default at start up is level 1. Cycle through the levels to get back to 1.
Pulse button (#1)	Short press	Triggers a light pulse.
Touch electrode (#9)	Short press	Triggers a light pulse when pressed together with the pulse button (#1).
Energy level indicator lights (green) (#5)	Constant	The number of lights indicates the energy level.
Ready/skin tone warning indicator light (green) (#7)	Constant	Device is ready for use.
	Blinks slowly	There is good contact with the skin.
	Blinks rapidly	Skin tone is too dark for safe use. No light pulses are emitted.
Bluetooth indicator light (blue) (#4)	Blinks	Bluetooth is active.

5.2 Skin colour sensor

Light-based hair removal on darker skin tones may result in adverse effects such as burns, blisters, and skin colour changes (hyper- or hypo-pigmentation). Silk'n Infinity has a skin colour sensor that detects the colour of your skin automatically. The skin colour sensor checks your skin colour at the start of, and occasionally during, each treatment.


- i** If you see no light pulse and the ready/skin tone warning indicator light (#7) blinks rapidly, this is an indication that your skin tone, measured by the skin colour sensor, is too dark for safe application. Try using the device on a different part of your body.

5.3 Energy levels

The energy level determines the intensity of the light pulse delivered to your skin. As the energy level increases, so do hair removal results, as well as the risk of possible side effects and complications.

- The device has 5 levels of energy from 1 (the lowest) to 5 (the highest), represented by 5 energy level indicator lights (#5).
- When the device is turned on, the energy level is set to level 1 automatically. Only one energy level indicator light (#5) will be on.
- To change the energy level, press the control button (#6) repeatedly to increase the energy level from 1 to 5. Once at level 5, press the control button (#6) again to reduce the energy level from 5 back to 1.

5.4 Pulsing and gliding

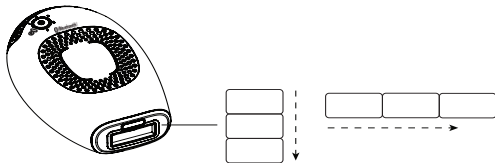Pulsing is the best technique for treating specific spots or small areas.

- Press the treatment surface (#2) against your skin. Work in rows (either across or down as shown in the image). Start at one end of each row and progress to the other end.
- The treatment surface (#2) creates temporary pressure marks on the treated area so you can see areas you have treated.

- After each pulse, move the treatment surface (#2) to the next spot. Keep doing this pulsing action for the area you are treating.
- Avoid gaps or overlaps as you use the device.

Gliding is the best technique for treating larger parts of the body.

Press the treatment surface (#2) firmly against your skin and slowly glide over your skin (either across or down as shown in the image). Make sure that the treatment surface (#2) is always touching your skin.

- At lower energy levels, the device takes only 1-1.5 seconds to recharge between pulses. This generates a nonstop sequence of flashes.
- Glide at a speed so that you avoid gaps or overlaps.
- ⚠ Do not treat the same area of skin more than once per hair removal session.
- ⚠ Stop use immediately if you have any pain.

5.5 Exposure to the sun

- The lowest energy level (1) is designed for safe use even after sun exposure.
- If you use energy levels 2-5, avoid further exposure to the sun for **4 weeks before and 2 weeks after** your treatment. This is because higher quantities of melanin expose the skin to greater risk when using any method of light-based hair removal. This applies to all skin types, even those that don't seem to tan quickly.
- We don't advise it, but if you have to go out in the sun after a treatment, be sure to apply sunscreen SPF30 or higher to the treated area. Do this for **2 weeks after** a treatment.

5.6 Exposure of pulsed light to the eyes

- Although you can use the device to treat facial hair, take extra care to avoid contact with the eyes, and use only below the cheekbones.
- 📌 The device will only emit repeated pulses when it is in firm contact with your skin, and you keep moving it across your skin.

6. What to expect with Silk'n Infinity

During a Silk'n Infinity treatment session it is normal to experience and feel:

- **A fan noise:** The cooling fan makes a noise similar to that of a hair dryer.
 - **A pop sound with each pulse:** When a pulse of eHPL™ light is activated, it is normal to hear a subtle pop sound simultaneously with the flash of light.
 - **A sensation of warmth and tingling:** During each galvanic pulse sequence it is normal to feel a mild sensation of warmth and tingling from the galvanic and light energies.
 - **Some mild red or pink colour:** During and just after your treatment session it is not uncommon to see some very mild, pink-like colour of the skin. This is usually most noticeable around the hairs themselves.
- ⚠ Stop use immediately if you see full redness of the skin, blistering or burns.

7. Possible side effects

When used according to the instructions, side effects and complications associated with the use of Silk'n Infinity are uncommon. However, every cosmetic procedure, including those designed for home use, involves some degree of risk. Therefore, it is important that you understand and accept the risks and complications that can occur with pulsed light hair removal systems designed for home use.

Minor skin discomfort

Home pulsed light hair removal is generally well-tolerated. However, some users feel mild discomfort during use. This is best described as a mild stinging sensation on the treated skin areas. The stinging sensation may last during or after treatment.

This means that you should use the device at a lower energy level, and if that doesn't help, you should stop using the device.

Skin redness

Your skin may become red immediately after using the device or within 24 hours. Redness generally clears up within 24 hours.

See your doctor if redness does not go away within 2 to 3 days.

Make sure you use the energy level that is comfortable for you and suitable for your skin colour.

Increased sensitivity of the skin

The skin of the treated area is more sensitive so you may encounter dryness or flaking of the skin.

Pigment changes

The device targets the hair shaft, in particular the pigmented cells in the hair follicle and the hair follicle itself. Nevertheless there is a risk of temporary hyper-pigmentation (increased pigment or brown discolouration) or hypo-pigmentation (whitening) to the surrounding skin. Risk of changes in skin pigmentation is higher for people with darker skin tones or heavily tanned skin. Usually discolouration or changes to skin pigment are temporary, and permanent hyper-pigmentation or hypo-pigmentation are rare.

Excessive redness and swelling

In rare cases, treated skin may become very red and swollen. This is more common in sensitive areas of the body. The redness and swelling should subside within 2 to 7 days and should be treated with frequent applications of ice. Gentle cleansing is OK, but avoid exposure to the sun.

Blisters or burns

Cool the affected area and apply anti-burn cream. Wait until the area has healed completely before another treatment, and use a lower energy level. **If the skin does not return to normal within 3 days, consult your physician.**

-  Stop use immediately if any of these effects are more severe.
-  See your physician if the effects do not go away within 2 to 3 days.

8. Treatment plan

A typical full hair growth cycle may take 18-24 months. During this time multiple Silk'n Infinity treatment sessions may be required in order to achieve permanent hair removal.

The efficiency of hair removal varies from person to person according to body area, hair colour and how the device is used.

Typical treatment plan during a full hair growth cycle:

-  **Treatments 1-4:** Plan two weeks apart.
-  **Treatments 5-7:** Plan four weeks apart.
-  **Treatments 8 +:** Treat as needed, until desired results are achieved.

Expected treatment outcome for hair removal

Fitzpatrick skin type ¹	Hair colour	Light fluence [J/cm ²]	Area of the body ²	Average number of treatments ¹	% Hair reduction
I-IV	Dark blond to brown	2.5-4	Lower or upper legs	Session 10-12, approx. four weeks apart.	60%
			Arms		50%
			Bikini Line		60%
			Underarms		55%
I-IV	Brown to black	2.5-4	Lower or upper legs	Session 8-10, approx. four weeks apart.	70%
			Arms		60%
			Bikini Line		70%
			Underarms		65%

1. Individual response depends on hair type as well as biological factors that may affect hair growth patterns. Some users may respond quicker or slower than the average number of treatments.

2. One cannot expect permanent hair removal in a single or even double treatment sessions. The duration of the resting period for hair follicles also depends on the body area.

9. Get started!

9.1 Connect the device

1. Remove the device and other components from the box.
2. Place the adapter into the socket (#10) and connect to the electricity supply. All the energy level indicator lights (#5) will go on and off repeatedly in a circular motion.

9.2 Do a spot test

It's important to do a spot test on each part of the body you wish to treat. A spot test checks that the skin colour is suitable for treatment, and that you know which energy level is comfortable for the area you are treating. Both skin colour and energy level may vary depending on the target treatment area.

1. To do a spot test, follow steps 1-5 and 7-9 in the next section.
2. Continue testing at increased energy levels until you have reached the highest level that you feel comfortable with. If within one hour you experience no adverse effects, you can start your full treatment at that energy level.
3. Repeat the spot test for each area you wish to treat.

9.3 Treatment steps

1. Remove all jewellery before use.
2. Before treatment, the skin should be shaved, clean, dry and free of any powders, antiperspirants or deodorants. **Do not wax, pluck or tweeze the hair out.**
3. Be sure to cover up tattoos, moles, warts or other marks on the skin that should not be treated. Use a material that does not absorb the light such as a white cloth or white medical tape.
4. Press the control button (#6). The device will turn on, the level 1 energy level indicator light (#5) will be on, and a fan sound will start.
5. Approximately 1 second after pressing the control button, the ready/skin tone warning indicator light (#7) (green) will turn on. The device is now ready for you to trigger the first pulse at the lowest energy level. Always start at energy level 1 for your first treatment.
6. Select the energy level you are most comfortable with.

- Place the treatment surface (#2) on the skin and press the pulse button (#1) to deliver one pulse. Make sure your hand touches the base electrode together with the pulse button (#1).

 If you do not make hand contact with the base electrode (#8) no galvanic treatment circulation will be emitted.

- Press the pulse button (#1).
 - The device first determines the skin colour. If the skin colour is light enough for safe application, the device will activate the galvanic circulation. You might feel a slight tingling sensation. This is completely normal. The galvanic treatment works to enlarge your pores for a better application of the treatment.
 - Directly afterwards, the device will flash a pulse of light onto your skin, and the ready/skin tone warning indicator light (#7) will turn off. You will see a bright flash of light and simultaneously hear a popping sound, which is a normal noise for the device. You will feel a mild sensation of warmth and tingling. The device will immediately recharge for the next pulse.
 - After a few seconds, when the device is ready to flash again, the ready/skin tone warning indicator light (#7) will turn on. The indicator light will start blinking when pressed against the skin.
- Remove the treatment surface (#2) from the treated area of skin.
- Decide if you wish to treat the area by **Gliding** or **Pulsing** (see section 5.4 Pulsing and gliding).
 - If you are **Pulsing**, place the treatment surface (#2) on an adjacent part of your skin, and trigger a pulse.
 - If you are **Gliding**, hold the pulse button (#1) down continuously and slide the device smoothly along the treated area.
- Keep **Pulsing** or **Gliding** until you have covered the target area.

 Try to avoid gaps or overlaps.

 If your skin blisters or burns, stop use immediately.

- When you have completed your treatment, long press the control button (#6), to switch off the device. Remember the last energy setting you used because the device always starts up at energy level 1.
- Disconnect the adapter from the electricity supply, then clean the device and store safely (see section 10. Cleaning & maintenance).

14. After your treatment session we recommend you apply a body lotion to the treated area to help your skin recover from the treatment faster.

10. Cleaning & maintenance

- Unplug the device.
 - To clean the device including the treatment surface (#2), wipe with a clean, dry cloth. Never use acetone or other abrasive fluids, because they will damage your device.
 - After cleaning, we recommend you store the device in its original box and keep it away from water.
- ⚠ Never immerse the device or any of its parts in water.

11. Connect your Silk'n Infinity to Bluetooth

Download and connect to Silk'n hair removal app.

The free Silk'n Hair Removal application supports your every step with Silk'n hair removal devices, creates your personalized treatment plan, helps you make the safest and most informed choices, and makes sure you follow through the entire course of treatment for best results.

As a Silk'n Infinity owner you can also enjoy Bluetooth connectivity. Just download the application and follow the registration process and Bluetooth pairing steps.

For a free download of the Silk'n hair removal app, scan:

 **Bluetooth**

12. Troubleshooting

Only authorized personnel are permitted to perform repairs. Do not attempt to modify, open or repair the device. You may expose yourself to dangerous electrical components and may cause serious harm to yourself. It will also void the warranty.

Problem	Checks
My device does not start.	<p>Check that:</p> <ul style="list-style-type: none"> the adapter is properly connected to the device and that it is plugged into the electricity supply.
There is no galvanic pulse when I press the pulse button (#1).	<p>Check that your hand is touching the base electrode when you press the pulse button (#1).</p> <p>The galvanic treatment circulation requires closing the galvanic circuit by pressing the pulse button (#1) and the electrode (#9) at the same time, in order to effectively treat your skin.</p>
There is no light pulse when I press the pulse button (#1).	<p>Check that:</p> <ul style="list-style-type: none"> the green ready/skin tone warning indicator light (#7) is blinking slowly. This indicates that you have good contact with the skin. the treatment surface (#2) is evenly and firmly pressed to the skin. For your safety, the pulse button (#1) will activate a pulse only if the treatment surface (#2) is firmly pressed against the skin. <p>If the ready/skin tone warning indicator light (#7) blinks rapidly, this is an indication that your skin tone, measured by the skin colour sensor, is too dark for safe use. Try using the device on a different body part or contact Silk'n customer support.</p>
One or more lights are on constantly.	<p>The device has malfunctioned. Unplug the device and do a restart. If the problem persists, contact Customer Service.</p>
The device has stopped working and is in standby mode.	<p>Press the control button (#6) again to start the device.</p>

13. Customer service

For more information see our website www.silk'n.eu. This manual is also available as a PDF download from the website. Contact Silk'n Customer Service to report unexpected operations or events, or any other problem with your device. Our contact details are:

Country	Service number	Email
Belgium	0800-29316	serviceconsommateurbe@silk'n.eu
France	0891-655557	serviceconsommateurfr@silk'n.eu
Germany	089 51 23 44 23	kundenservicede@silk'n.eu
Italy	+31(0)180-330550	servizioclientiit@silk'n.eu
Netherlands	0900-2502217	klantenservicenl@silk'n.eu
Spain	900 823 302	servicioalcliente@silk'n.eu
UK	0906-2130009	customerareuk@silk'n.eu
Other countries	+31 (0)180-330 550	info@silk'n.eu

14. Warranty

Refer to your Silk'n warranty card.

15. Disposal

Sort the packaging before you dispose of it. Dispose of paperboard and cardboard with the recycled paper service and wrappings with the appropriate collection service.

(Applicable in the European Union and other European countries with separate collection systems of recyclable materials). WEEE-Reg.Nr. DE 31619958. **Old appliances may not be disposed of in household waste!** Should the product no longer be capable of being used at some point in time, **dispose of it in accordance with the regulations in force in your city or county.** This ensures that old appliances are recycled in a professional manner and also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

*Products labelled with Cd = cadmium, Hg = mercury, Pb = lead.

16. Specifications

Model	H3103
Galvanic	70-400 Microampere (μA)
Spot size	0.9cm x 3cm [2.7cm ²]
Technology	Enhanced Home Pulsed Light™
Max energy level	Max 4J/cm ²
Wavelength	520-1200nm
Pulse light duration	500-800 μs
Mode of operation	Continuous
Operation and safety	Skin tone sensor ensures use only on appropriate skin types
Repetition rate	1.0Hz maximum
System weight	225gr
Transport & storage between uses and storage condition	Temperature: -40° to +70°C Relative humidity: 10% - 90%rH Atmospheric pressure: 500 - 1060hPa
Operating condition	Temperature: 10° - 30°C Relative humidity: 30% - 75%rH Atmospheric pressure: 700 - 1060hPa
Time from minimum storage temperature between uses to minimum operating temperature with ambient temperature of 20°C	30 minutes
Time from maximum storage temperature between uses to maximum operating temperature with ambient temperature of 20°C	30 minutes
Protection against electric shock	Class II medical equipment type BF applied part
Adapter	KSA-60B-200300D5
	Input: 100-240V; 50-60Hz; 1.8A Output: 20Vdc; 3A

Pulse parameters (seconds)	Energy level 1: 0.5 Energy level 2: 0.7 Energy level 3: 0.9 Energy level 4: 1.0 Energy level 5: 1.3
Bluetooth LE Operating Frequency Transmitter Power	2,402 MHz ~ 2,480 MHz Below 3mW

17. Legal notice

Copyright 2023 © Silk'n Beauty Ltd. All rights reserved.

Home Skinovations Ltd. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd. is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd. No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification. Silk'n Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silk'n Beauty Ltd. Specifications are subject to change without notice. Silk'n and the Silk'n logo are registered trademarks of Silk'n Beauty Ltd., Invention Works B.V. Galileistraat 17, 3029 AL Rotterdam, The Netherlands www.silkn.eu, info@silkn.eu.

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Copyright 2023 © Silk'n Beauty Ltd. Alle Rechte vorbehalten.
Copyright 2023 © Silk'n Beauty Ltd. Tous droits réservés.
Copyright 2023 © Silk'n Beauty Ltd. Alle rechten voorbehouden.
Copyright 2023 © Silk'n Beauty Ltd. Todos los derechos reservados.
Copyright 2023 © Silk'n Beauty Ltd. Todos os direitos reservados.
Copyright 2023 © Silk'n Beauty Ltd. Tutti i diritti riservati.